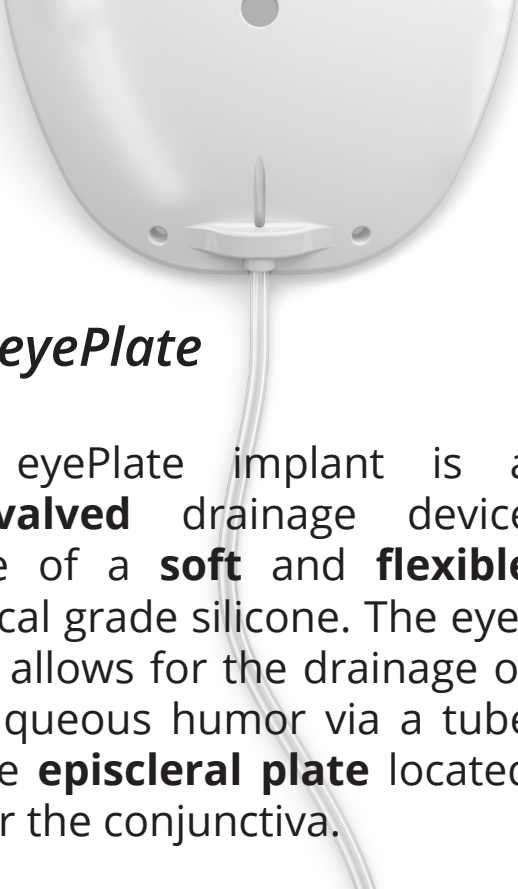




The eyePlate

The eyePlate implant is a **non-valved** drainage device made of a **soft** and **flexible** medical grade silicone. The eye-Plate allows for the drainage of the aqueous humor via a tube to the **episcleral plate** located under the conjunctiva.

Thanks to its unique design **avoiding muscle isolation**, the surgical procedure can be done under **local anesthesia**.

The eyePlate can be used either as a **stand-alone device** or in conjunction with the eyeWatch device.

The eyePlate-200
& eyePlate-300
are CE marked products



RHEON
MEDICAL 

Revolutionizing the
surgical treatment of
glaucoma.

Order now!

website: www.rheonmedical.com
email: info@rheonmedical.com

Rheon Medical
EPFL Innovation Park-B
1015 Lausanne
Switzerland



eyePlate
Glaucoma Drainage Device

RHEON
MEDICAL 

www.rheonmedical.com

Information leaflet

The eyePlate

**Soft. Flexible.
Ideal design.**

Tubing

With a lumen of 0.3 mm and external diameter of 0.63 mm

Smooth and thin

Smooth and thin (0.8mm) medical grade silicone, increasing its flexibility. Can be folded and placed under the conjunctiva.

Exists in 2 sizes

200mm² & 300mm²

Optimized design

The endplate has been designed to maximize surface area while keeping the shape narrow, to limit the contact with the recti muscles.

